

PRODUCT OVERVIEW

This electric heating pad provides adjustable surface heat for temporary comfort and relaxation of tired, stiff, or sore muscles. Use only as directed and never use while sleeping.

SPECIFICATIONS

- Rated supply: 120 V AC, 60 Hz
- Rated power: 50 W
- Heat settings: Low, Medium, High
- Automatic shut-off: 60 minutes
- Soft-touch removable, washable cover

PACKAGE CONTENTS

- Heating pad
- Power cord and controller
- User guide

PRODUCT INFORMATION

Keep this guide for future reference. Read all instructions before using the product.

CUSTOMER SUPPORT

Distributed by United Yoram Inc.
Brooklyn, NY 11236
drfranklyns.com
Instagram: @drfranklyns

DOCUMENT CONTROL

Brand: Dr. Franklyn's | Model: HP00XF | Revision: 1.0 | Issue date: September 2026

SCAN FOR PRODUCT DETAILS



drfranklyns.com


IMPORTANT SAFETY INSTRUCTIONS

- Do not use while sleeping.
- Do not use on infants, children, incapacitated persons, or anyone insensitive to heat.
- Do not sit or lie on the pad. Do not fold, bunch, crease, or use while compressed.
- Keep dry. Do not use in wet or moist environments.
- Do not pierce with pins, needles, or sharp objects.
- Inspect the pad, controller, connector, and cord before use. Do not use if damaged or malfunctioning.
- Do not trap the cord, wrap it tightly, or allow it to contact hot surfaces.
- Unplug when not in use and before cleaning.
- Stop use if the heat feels excessive or causes burning, pain, redness, or skin irritation.
- Consult a healthcare professional before use if you have diabetes, poor circulation, reduced sensation, or a serious medical condition.

For household use only. Use only the supplied controller and connectors.



HEALTH | SUPPORT | RECOVERY

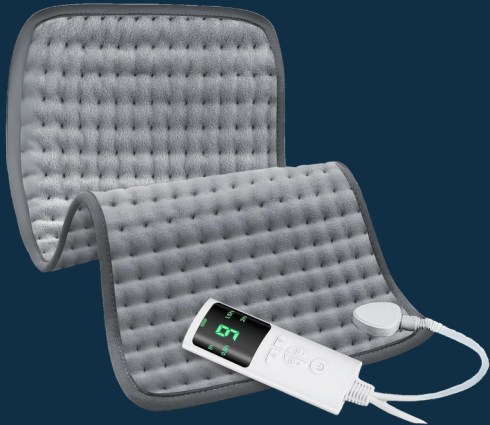

MODEL HP00XF

HEATING PAD

Adjustable heat for everyday comfort

READ INSTRUCTIONS COMPLETELY BEFORE USE

SET UP AND USE


- 1 UNFOLD**
Fully unfold the pad and place it flat over the intended body area. Never place your body weight on the pad.
- 2 CONNECT**
Align and firmly connect the controller connector to the heating-pad connector before plugging into the outlet.
- 3 POWER**
Plug into a properly installed 120 V AC outlet and press the power button.
- 4 SELECT**
Choose Low, Medium, or High. Start at the lowest comfortable setting and monitor the skin frequently.
- 5 FINISH**
The pad shuts off automatically after 60 minutes. Unplug after use and allow it to cool before storing.

INDICATOR GUIDE


- Low: green indicator
- Medium: yellow indicator
- High: red indicator

Check the skin regularly. Stop use immediately if the pad feels too hot or causes discomfort, persistent redness, or irritation.

CONTROLLER AND CONNECTIONS

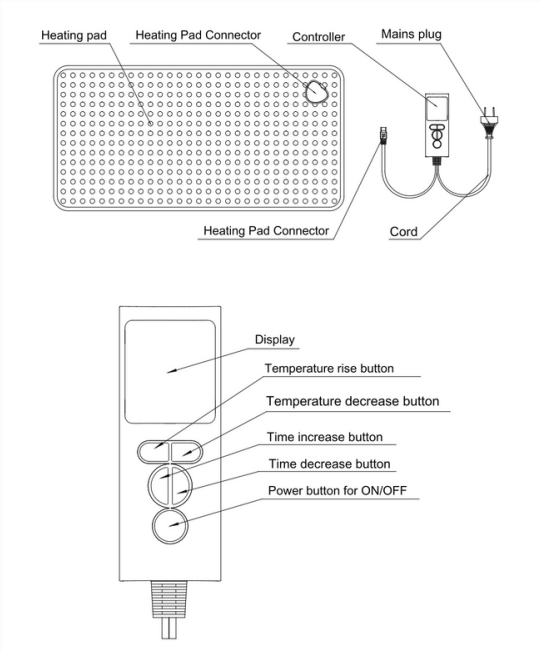

CONNECT BEFORE POWERING

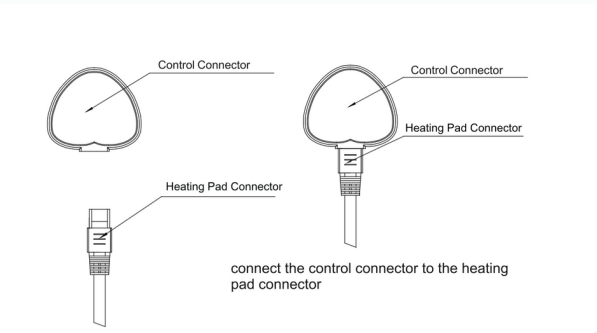

CLEANING

- Unplug and disconnect the controller before cleaning.
- Remove the outer cover and hand wash or machine wash it in cold water.
- Lay the cover flat to air dry completely. Do not machine dry.
- Do not wash or immerse the internal heating element, controller, cord, or connectors.
- Reconnect only after every part is completely dry.

STORAGE

- Unplug and allow the pad to cool completely.
- Store flat in a clean, dry place.
- Do not sharply fold the pad or wrap the cord tightly around the controller.
- Do not place heavy objects on the pad during storage.

TROUBLESHOOTING

NO POWER

Confirm the pad and controller connectors are fully seated, the outlet works, and the power button has been pressed.

NO HEAT

Select a heat setting and allow several minutes to warm. If it remains cold, unplug and discontinue use.

TURNS OFF

The 60-minute automatic shut-off is normal. Allow the unit to cool before restarting if continued use is appropriate.

DAMAGE OR ERROR

Unplug immediately. Do not attempt to repair the pad, controller, cord, or connector.