

QUICK REFERENCE

Brand / model	Dr. Franklyn's ITW-1
Power source	Integrated rechargeable battery
Pad size	400 × 200 × 4.95 mm
LED array	120 × 0.5 W LEDs
Wavelengths	660 nm red + 850 nm NIR
Irradiance	Over 100 mW/cm² at 0 in
Automatic timer	20 minutes
Intensity	50% / 75% / 100%
Operating modes	Steady / breathe / 10 Hz / 40 Hz
Charging	Type-C
Net weight	Approximately 0.58 kg
Material	Flexible fabric
LED life	Over 50,000 hours

CORDLESS OPERATION

Charge the integrated battery fully before first use. Once charged, disconnect the charging cable and use the pad cordlessly. Recharge only through the device charging port with the supplied charging cable.



CHARGE



DISCONNECT




USE CORDLESS

CUSTOMER SUPPORT

Distributed by United Yoram Inc.
Brooklyn, NY 11236 • drfranklyns.com
Instagram: @drfranklyns

SCAN FOR PRODUCT DETAILS

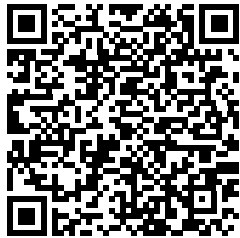

drfranklyns.com

READ BEFORE USE

This flexible red and near-infrared light pad is intended for personal wellness use on intact skin. Read all instructions before use and keep this guide for future reference.

IMPORTANT

Do not stare at illuminated LEDs. Stop use if you experience pain, burning, dizziness, headache, unusual skin changes, or worsening symptoms.

DO NOT USE

- If pregnant, unless cleared by a healthcare professional.
- If you have epilepsy or a diagnosed light-sensitive condition.
- If you take medication that may increase photosensitivity, unless cleared by a qualified professional.
- Over open wounds, infected skin, active bleeding, or areas with impaired sensation.
- Near water, while bathing, or when the pad, charging port, or charging cable is wet.

CARE AND STORAGE

- Switch off, disconnect the charging cable, and allow the pad to cool before cleaning.
- Wipe gently with a soft, lightly damp cloth. Do not use harsh chemicals.
- Do not immerse, machine wash, tumble dry, iron, fold sharply, or place heavy objects on the pad.
- Store flat or loosely rolled, away from direct sunlight, moisture, excessive heat, children, and pets.

WHEN TO SEEK ADVICE

If discomfort is sudden, severe, unexplained, or does not improve, stop using the device and consult a healthcare professional. This product does not replace medical diagnosis or treatment.



RED LIGHT THERAPY PAD

RED + NEAR-INFRARED WELLNESS



MODEL ITW-1

FLEXIBLE • PORTABLE • ADJUSTABLE

USER MANUAL

660 nm red light • 850 nm near-infrared light
Read all safety instructions before use

1. ABOUT THE THERAPY PAD

The ITW-1 combines 660 nm visible red light with 850 nm near-infrared light. Red light is used for surface-level wellness applications, while near-infrared light reaches deeper tissue. Individual results vary.


2. WHAT'S INCLUDED

- Cordless red light therapy pad
- Adjustable elastic strap
- Charging cable
- User manual



3. COMMON APPLICATION AREAS

The flexible pad can be positioned around the back, shoulders, knees, feet, arms, and other suitable external areas. Place the illuminated surface toward the treatment area.

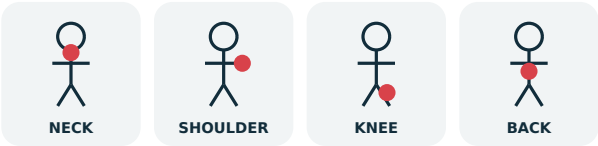

COMFORTABLE PLACEMENT

Secure snugly enough to maintain position, but never so tightly that circulation is restricted.

4. HOW TO USE

- 1

INSPECT
Check the pad, charging port, and charging cable. Do not use damaged or wet components.
- 2

CHARGE
Before cordless use, fully charge the integrated battery with the supplied charging cable. Disconnect the cable when charging is complete.
- 3

POSITION
Place the LED side over the desired area and secure with the adjustable strap.
- 4

POWER ON
Press the power button. The device display or indicator will illuminate.
- 5

SELECT SETTINGS
Choose the operating mode and intensity level using the integrated control panel.
- 6

COMPLETE SESSION
The device switches off automatically after 20 minutes. Switch off after use and recharge when needed.

20-MINUTE AUTOMATIC SESSION


Each operating cycle is timed for 20 minutes. The pad switches off automatically when the cycle ends. Stop earlier if the pad feels uncomfortable or too warm.



CHARGING NOTE

Do not use the pad while it is charging. Keep the charging port and cable dry and free of damage.

5. INTEGRATED CONTROL PANEL

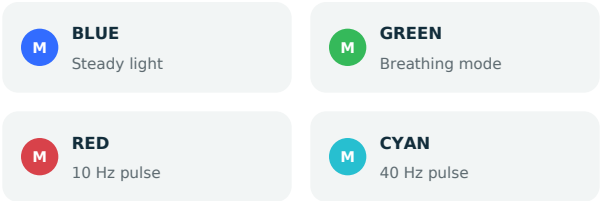

Timer	Fixed 20-minute session with automatic shutoff.
Charging	Type-C charging port.
Battery	Red indicator: low battery. Green: full / working.
Intensity	Blue: 50%. Green: 75%. Red: 100%.

6. TROUBLESHOOTING

- No power or display**
Charge the integrated battery, disconnect the charging cable, and press the power button again.
- LEDs do not illuminate**
Confirm the selected mode and intensity. Switch off, wait 10 seconds, and restart.
- Pad feels too warm**
Switch off, remove the pad, and allow it to cool. Use a shorter session and inspect the skin.
- Controls do not respond**
Switch off for 30 seconds, then restart. Recharge if needed. Stop use if the problem continues.

7. DISPOSAL

Dispose of electrical components and accessories according to local electronic-waste requirements. Do not burn or place damaged electrical components in household trash where prohibited.