



CUSTOMER SUPPORT

Distributed by United Yoram Inc. | Brooklyn, NY 11236
drfranklyns.com
support@drfranklyns.com
Instagram: @drfranklyns

PRODUCT INFORMATION

Product: Dr. Franklyn's Pneumatic Compression Therapy Device
Model / SKU: DVT-3
Manual revision: 1.1 | September 2026
Limited warranty: 1 year
Expected product lifespan: 3 years

KEEP THIS GUIDE

Keep this manual for future reference. Follow the prescribing clinician's directions whenever they differ from general guidance in this booklet.

IMPORTANT SAFETY INFORMATION

The Dr. Franklyn's Pneumatic Compression Therapy Device applies intermittent external pressure to support venous blood flow. When prescribed by a qualified clinician, it may be used as part of a plan to help reduce the risk of deep vein thrombosis (DVT). It is not intended to treat an existing blood clot.

- Do not use if deep vein thrombosis, thrombophlebitis, pulmonary embolism, severe arterial disease, acute infection, or uncontrolled heart failure is known or suspected unless specifically directed by a qualified clinician.
- Do not apply over open wounds, recent skin grafts, active cellulitis, gangrene, severe edema, or severe deformity without medical approval.
- Begin at the lowest comfortable pressure. Never tighten the wrap enough to restrict circulation or cause pain, numbness, tingling, or skin-color changes.
- Keep the controller, USB cable, battery compartment, and connections dry. Do not immerse or operate near water.
- Do not disassemble, modify, drop, crush, or expose the controller or batteries to excessive heat, flame, or flammable anesthetic mixtures.

DO NOT USE / ASK A CLINICIAN

- Ask a clinician before use with neuropathy, impaired sensation, significant vascular impairment, recent surgery, pregnancy, cancer, or an implanted medical device.
- Children and people unable to communicate discomfort require direct caregiver and clinician supervision.

STOP USE

Stop immediately and seek medical advice for chest pain, shortness of breath, sudden limb swelling, increasing pain, numbness, skin discoloration, dizziness, or any unusual reaction.



PNEUMATIC COMPRESSION THERAPY DEVICE

INTERMITTENT AIR COMPRESSION FOR CIRCULATION SUPPORT



MODEL DVT-3

USER MANUAL

READ ALL INSTRUCTIONS BEFORE USE

HEALTH | SUPPORT | RECOVERY

INTENDED USE AND MEDICAL PURPOSE

The Dr. Franklyn's Pneumatic Compression Therapy Device is intended to provide intermittent external compression to help support venous circulation, reduce temporary limb fatigue and swelling, and promote recovery. Under clinician direction, it may be used as part of DVT-prevention care for appropriate patients.

IMPORTANT

This device does not diagnose or treat an existing blood clot. Suspected DVT or pulmonary embolism requires immediate medical evaluation.

HOW IT WORKS

The controller cyclically inflates and releases the connected wrap at the selected pressure and program. The repeated pressure supports movement of venous blood through the treated limb.

KEY FEATURES

- Five selectable pressure levels from 20 to 100 mmHg.
- Two air ports for one or two compatible wraps.
- Selectable intermittent compression programs.
- Universal wrap for the leg, calf, thigh, arm, elbow, shoulder, or waist/lower-back area.



SETUP AND OPERATION

- 1. Inspect the controller, wrap, air bladder, tubing, ports, straps, and battery cover. Do not use damaged components.
- 2. Position the wrap over the selected area and fasten it snugly without restricting circulation.
- 3. Connect the tubing firmly to air port P1. P2 may support a second compatible wrap.
- 4. Press and hold the Power button to turn on the controller.
- 5. Select a comfortable pressure from 20 to 100 mmHg. Begin at the lowest setting.
- 6. Press the Mode button to select the desired intermittent compression program.
- 7. Monitor skin color, warmth, sensation, and comfort throughout the session. When finished, turn off, loosen the wrap, and disconnect the tubing.

QUICK REFERENCE


Item	Specification
Pressure range	20-100 mmHg
Pressure levels	5
Air ports	2 (P1 and P2)
Power	Two protected flat-top 18650 Li-ion batteries; not included
Battery capacity	1800-2600 mAh each
Battery size	65-68.5 mm long x 18-18.6 mm diameter

PRESSURE AND PORT GUIDANCE

- Increase pressure only as needed and directed. Compression should feel firm and rhythmic, never painful.
- Do not kink, pinch, or sharply bend the air tubing. Stop if inflation becomes uneven or uncomfortable.

BATTERY INSTALLATION

Requires two matching protected flat-top 18650 rechargeable lithium-ion batteries. Batteries are sold separately and are not included.

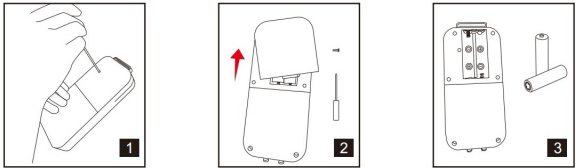

- 1. Disconnect the USB cable. Use the supplied screwdriver to loosen and remove the battery-cover screw.
- 2. Slide the upper battery cover upward and remove it.
- 3. Insert both batteries according to the + and - polarity markings.
- 4. Reinstall the cover and securely tighten the screw before use.

CARE AND STORAGE

- Turn off and disconnect the controller before cleaning.
- Wipe the controller with a soft, lightly damp cloth. Never immerse it.
- Wipe the wrap gently and allow it to air dry completely. Do not machine wash or tumble dry unless the product label specifically permits it.
- For long-term storage, remove the batteries. Store all components clean and dry, away from sunlight, heat, moisture, children, pets, and sharp objects.

TROUBLESHOOTING

Issue	Action
Will not turn on	Check battery polarity and charge; secure the battery cover.
Weak or uneven pressure	Refit the wrap and check tubing, ports, and connections.
Abnormal operation	Stop use, power off, and contact customer support.

WHAT IS INCLUDED

DVT-3 controller, universal pneumatic wrap, air tubing, USB charging cable, screwdriver, and this user manual. Batteries are not included.

