

## IMPORTANT SAFETY

- Keep away from open flame, cigarettes, heaters, and other heat sources. Polyurethane foam can burn.
- Not intended for infants or anyone unable to reposition independently without appropriate supervision.
- Make sure the pad lies flat and does not create a trip or transfer hazard at the mattress edge.
- Discontinue use if the foam tears, crumbles, becomes wet, develops odor, or causes irritation.

## BEFORE EACH USE

- Check that the pad is clean, dry, flat, and centered.
- Inspect for tears, compression, discoloration, or breakdown.
- Confirm that fitted bedding does not bunch or pull the pad out of position.

## IMPORTANT NOTE

This comfort pad is not a replacement for a prescribed pressure-redistribution mattress or a clinician-directed pressure-injury prevention plan.



## PRODUCT INFORMATION

|               |                             |
|---------------|-----------------------------|
| PRODUCT       | Egg Crate Foam Mattress Pad |
| SIZE          | Twin                        |
| SKU           | MAT-1                       |
| MATERIAL      | Polyurethane foam           |
| WARRANTY      | 1 year                      |
| EXPECTED LIFE | 3 years                     |

## INCLUDED

1 Twin Egg Crate Foam Mattress Pad

## SCAN FOR PRODUCT DETAILS



Product information, support, and updates.


[drfranklyns.com](https://drfranklyns.com)

United Yoram Inc. | Brooklyn, NY 11236 | [support@drfranklyns.com](mailto:support@drfranklyns.com)



# Egg Crate Foam Mattress Pad

TWIN | CUSHIONING AND AIRFLOW

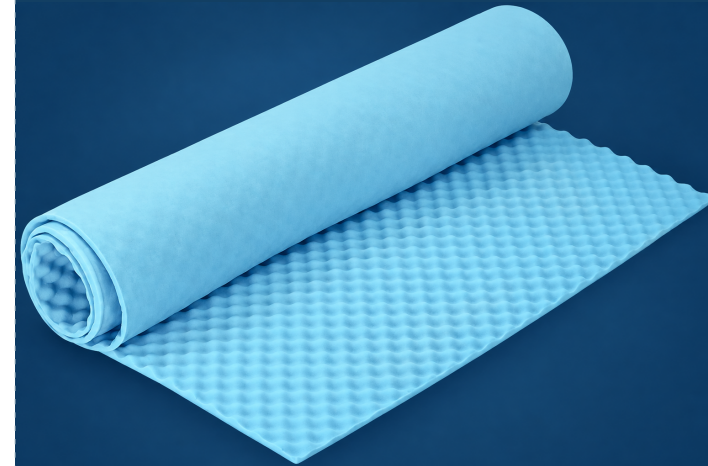

USER MANUAL | MAT-1

READ ALL INSTRUCTIONS BEFORE USE

HEALTH | SUPPORT | RECOVERY

## PRODUCT OVERVIEW

The Dr. Franklyn's Egg Crate Foam Mattress Pad adds a lightweight cushioning layer to a twin mattress. Its textured surface is designed to help distribute body weight, improve airflow, and reduce localized pressure for everyday comfort.



## KEY FEATURES

- Egg-crate surface creates multiple cushioning zones.
- Open texture encourages airflow beneath the body.
- Adds comfort to firm or aging mattresses.
- Lightweight design is easy to position and remove.
- Provides a protective layer that may help reduce mattress wear.

## SETUP AND USE

1

### UNROLL

Remove all packaging carefully. Do not use a knife near the foam.

2

### ALLOW TO EXPAND

Lay the pad flat in a well-ventilated area and allow it to regain its natural shape before use.

3

### ORIENT

Place the textured egg-crate peaks facing upward unless your healthcare professional directs otherwise.

4

### CENTER

Position the pad evenly on top of the twin mattress. Keep all edges flat and within the mattress perimeter.

5

### COVER

Place a clean fitted sheet over the pad and mattress. Avoid very tight bedding that compresses the foam.

6

### CHECK

Before lying down, confirm the pad has not shifted, folded, or bunched.

### POSITIONING TIP

Re-center the pad whenever bedding is changed. For users at risk of skin breakdown, follow the turning, repositioning, and skin-check plan provided by a clinician.

## CLEANING AND CARE

- Remove bedding and gently vacuum loose debris using low suction.
- Spot clean only with a soft cloth, mild soap, and cool water. Do not saturate the foam.
- Press with a dry towel to remove moisture, then air dry completely in a ventilated area.
- Do not machine wash, dry clean, bleach, iron, or use a heated dryer.
- Do not replace bedding or use the pad until completely dry.

## STORAGE

- Store clean and completely dry.
- Keep flat or loosely rolled; do not sharply fold or compress for long periods.
- Protect from moisture, direct sunlight, heat, sharp objects, and heavy items.

## REPLACE THE PAD IF

- Foam remains flattened and no longer rebounds.
- Cracks, tears, crumbling, permanent staining, or persistent odor appear.
- The pad no longer stays flat or supportive.

## COMFORT AND SKIN CHECKS

Stop use if discomfort or skin irritation develops. Individuals with limited mobility, reduced sensation, or existing skin concerns should use this product only with appropriate clinical guidance and regular skin inspection.