




CUSTOMER SUPPORT

Distributed by United Yoram Inc. | Brooklyn, NY 11236
drfranklyns.com
support@drfranklyns.com
Instagram: @drfranklyns

PRODUCT INFORMATION

Product: Dr. Franklyn's Advanced Wrist Support
Model / SKU: WS-C
Manual revision: 1.1 | September 2026
Limited warranty: 1 year
Expected product lifespan: 3 years

SCAN FOR PRODUCT DETAILS

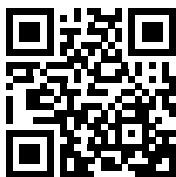

drfranklyns.com

KEEP THIS GUIDE

Keep this manual for future reference. Follow the prescribing clinician's directions whenever they differ from general guidance in this booklet.

IMPORTANT SAFETY INFORMATION

Use the Dr. Franklyn's Advanced Wrist Support only as described in this guide or as directed by a qualified healthcare professional.

- Do not over-tighten. The brace must not restrict circulation or cause numbness, tingling, coldness, swelling, or skin-color changes.
- Do not apply over open wounds, infected skin, or unprotected surgical incisions.
- Inspect the splint, straps, stitching, thumb opening, and hook-and-loop closures before every use. Do not use a damaged brace.
- Do not cut, heat, bend, or otherwise modify the support splint unless instructed by a qualified clinician.
- Remove the brace before bathing or swimming and keep it away from direct heat and open flame.

DO NOT USE / ASK A CLINICIAN

- Seek medical evaluation for suspected fracture or dislocation, severe swelling, deformity, loss of movement, or sudden severe pain.
- Ask a clinician before use after surgery or with poor circulation, impaired sensation, fragile skin, or a significant skin condition.

STOP USE

Remove the support and contact a healthcare professional if pain increases or if numbness, tingling, weakness, discoloration, swelling, pressure injury, or persistent skin irritation develops.



ADVANCED WRIST SUPPORT

STABILIZATION FOR DAILY COMFORT & RECOVERY



MODEL WS-C

USER MANUAL

READ ALL INSTRUCTIONS BEFORE USE

HEALTH | SUPPORT | RECOVERY

INTENDED USE

The Dr. Franklyn's Advanced Wrist Support is designed to stabilize and support the wrist while helping limit unwanted movement during daily activity, rest, and clinician-directed recovery. The universal design may be worn on the left or right wrist.

COMMON APPLICATIONS

- Mild to moderate wrist sprains and strains.
- Tendon irritation, repetitive strain, and overuse discomfort.
- Carpal-area support when recommended by a healthcare professional.
- Weak or unstable wrists and post-injury rest.

KEY FEATURES

- Reinforced splint support for structured stabilization.
- Three-strap adjustable system for a personalized fit.
- Breathable interior for comfortable daily wear.
- Thumb-hole design and quick hook-and-loop application.

WEAR SCHEDULE

Follow your clinician's recommended wear time. The brace should feel secure and supportive, never painfully tight.

HOW TO APPLY

Apply to clean, dry skin. Keep the wrist in a comfortable neutral position throughout application.

- 1. Fully open all three straps and identify the thumb opening and reinforced splint.
- 2. Slide the hand into the brace and guide the thumb through the opening.
- 3. Position the splint along the wrist so the brace is centered and does not press into the palm crease.
- 4. Fasten the wrist strap, then secure the remaining straps evenly toward the forearm.
- 5. Adjust each strap until the support is snug and stable without restricting circulation.
- 6. Move the fingers and thumb. Confirm that the brace does not pinch and the wrist remains comfortably supported.

FIT CHECK

Check	Correct fit
Position	Brace is centered; thumb opening lies flat.
Support	Wrist feels stable in a comfortable neutral position.
Circulation	Fingers remain warm with normal color and sensation.
Straps	Secure and flat, without painful pressure or bunching.

REMOVAL

- Release all straps gradually, then slide the hand out without twisting the wrist.
- Check the skin after removal. Persistent redness or irritation requires a fit review.

CLEANING AND CARE

- Close hook-and-loop fasteners before cleaning to reduce snagging.
- Remove the splint only if your supplied brace is designed for splint removal.
- Hand wash in cool water using mild soap. Rinse thoroughly and air dry flat.
- Do not machine wash, bleach, iron, tumble dry, dry clean, or place near direct heat.
- Wipe a removable splint with a soft damp cloth and allow it to dry completely before reinserting.

STORAGE

Store clean and dry with straps loosely closed. Keep away from direct sunlight, excessive heat, moisture, pets, and sharp objects.

WHAT IS INCLUDED

One Dr. Franklyn's Advanced Wrist Support with reinforced splint and three adjustable straps, plus this user manual.

DAILY CHECK

Inspect both the brace and your skin before and after each use. Replace the support if the splint, fabric, stitching, or closures become damaged.